

### **FUTURE BEACON**

https://future-beacon.co











## PRESS IV

#### (WITH PRE COMPRESSION)

**DOUBLE SIDED HIGH SPEED ROTARY TABLET PRESS** 

45 Stn. 'D', 61 Stn. 'B', 75 Stn. 'BB'

25 Stn. Slugging Type (For veterinary Tablets) (40mm dia) 35 Stn. Slugging Type (For veterinary Tablets) (30mm dia)





#### **SALIENT FEATURES**

- It is a high speed double rotary tablet press for higher tablet output.
- Machine is having a square design with full consideration of cGMP standards.
- The machine is having PRE-COMPRESSION facility.
- The machine is having A.C. variable frequency drive for the main motor.
- The machine is having Electromagnetic clutch arrangements.
- All the tablet parameters can be fine tuned during the operation of the machine from outside the tableting zone, with front controls at outside of the machine, having SEIKO type dials.
- The machine is having A.C. variable frequency drive for force feeding motors, which enables operators to synchronize machine's speed & feed both.
- The machine having Motorized Hydraulic Powerpack for pressure loading and overload release.
- The machine is having Helical gear ring & pinion gear, to drive the turret.
- All upper Guards of transparent acrylic material with interlock switches for safety.
- Lower Guards of polished Stainless Steel SS 304.
- The machine is having excellent BOTTOM DRIVEN FORCE FEEDING SYSTEM for material feeding.
- The machine is having three piece turret with SS-316 middle plate.
- Electro less nickel plating is provided to the turret.
- Two layer attachment can also be provided, (Optionally).
- PLC with color MMI touch screen (Optionally).





#### With Auto Force Rejection Through Compaction Force

- With load cell at Main Compression for Pressure measurement.
- With auto sampling based on, either time or revolution of punch no.
- With auto rejection of any single tablet out of prescribed compaction force limit.
   SCADA with industrial PC. having 15" colour touch screen (21 CFR part 11 compliance)
- With graphical as well as numeric representation of each individual punch force.













# PRESS IV 45 Stn. 'D', 61 Stn. 'B', 75 Stn. 'BB 25 Stn. Slugging Type (For veterinary Tablets) (40mm dia) 35 Stn. Slugging Type (For veterinary Tablets) (30mm dia)

#### (WITH PRE COMPRESSION)

**DOUBLE SIDED HIGH SPEED ROTARY TABLET PRESS** 


Machine Type	PRESS IV-45	PRESS IV - 61	PRESS IV - 75	PRESS IV - 25 - SLG	PRESS IV - 35 - SLG
Number of Stations	45	61	75	25	35
Type of Tooling	D	В	BB	SLUGGING	SLUGGING
Output-Tablets/Hr.	65,000 to	87,840 to	1,08,000 to	(Min.) 36,000 /	(Min.) 50,400 /
(Min./Max.)*	3,25,000	4,39,200	5,40,000	(Max.) 1,05,000	(Max.) 1,45,000
Max. Operating Pressure(Main)	10 Tons	6.5 Tons	6.5 Tons	12 Tons	12 Tons
Max. Operating Pressure(Pre)	1 Ton	1 Ton	1 Ton	01 Tons	01 Tons
Depth of fill	20 mm	17.4 mm	17.4 mm	30 mm	30 mm
Max. Tablet Dia	25 mm	16 mm	11.1 mm	40 mm	30 mm
Upper Punch enrty	1.5 to 8 mm	1.5 to 8 mm	1.5 to 8 mm	1.5 to 8 mm	1.5 to 8 mm
Power Total	8.75 kw/11.7 HP	8.75 kw/11.7 HP	8.75 kw/11.7 HP	12.48 Kw / 16.7 HP	12.48 Kw / 16.7 HP
Power Main Electric Motor	11.25 kw /15 HP	11.25 kw /15 HP	11.25 kw /15 HP	11.25 kw /15 HP	11.25 kw /15 HP
	1440 RPM, 415 V.	1440 RPM, 415 V.	1440 RPM, 415 V.	960 RPM, 415 V.	960 RPM, 415 V.
	50 Hz, 3 Ph.	50 Hz, 3 Ph.	50 Hz, 3 Ph.	50 Hz, 3 Ph.	50 Hz, 3 Ph.
Power	0.19 KW / 0.25 HP	0.19 KW 0.25 HP	0.19 KW / 0.25 HP	0.19 KW / 0.25 HP	0.19 KW / 0.25 HP
(Feeder Motors 2 Nos.)	1440 RPM, 415 V.	1440 RPM, 415 V.	1440 RPM, 415 V.	1440 RPM, 415 V.	1440 RPM, 415 V.
	50 Hz, 3 Ph.	50 Hz, 3 Ph.	50 Hz, 3 Ph.	50 Hz, 3 Ph.	50 Hz, 3 Ph.
Overall Dimension (mm)	Guard Open2500 (L) X 2500 (W) X 2275 (H)			2570 (L) x 2800 (W) x	2570 (L) x 2800 (W) x
				2330 (H)	2330 (H)
Overall Dimension (mm)	Guard Close1500 (L) X 1580 (W) X 1900 (H)			1500 (L) X 1580 (W) X	1500 (L) X 1580 (W) X
				1900 (H)	1900 (H)
Net Weight (Approx)	3800 kg			3800 Kg.	3800 Kg.
Case Dimension (mm)	3100 (L) X 3200 (W) X 2450 (H)			3100 (L) X 3200 (W) X	3100 (L) X 3200 (W) X
				2450 (H)	2450 (H)
Gross Weight (Approx)	4150 kg			4150 kg	4150 kg


Pre-pressure and Main pressure could be changed depending on the size of the compress and filling depth, powder or granulation. It is dependent on the physical specification of the used material for the compress.

The above specification is the subject to change without prior notice for the technical development.

It is Recomended not to use large depth of fill, minimum tablet thickness and maximum compression pressure for maximum

performance consistently. Technical rights to change catalogues, specifications and features are reserved with the manufacturer. Actual output may vary against indicated output on the product under compression and RPM on lower side.

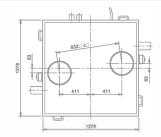

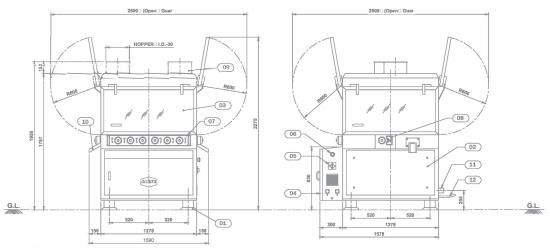





Address: Queen Rania St 337, Amman

Phone: +962(78)9990990

Email: info@futurebeacon.co

Website: https://future-beacon.co